

# Vico dressing table

**12207** The Vico dressing and examination table has been designed to dress or examine patients quickly and efficiently. Numerous possibilities and accessories make the work of the nursing staff easier, whereas the patient's welfare is increased.

## Carriage

- Hydraulic high/low adjustment by means of support arms.
- The steel frame is finished with a removable vacuum-formed synthetic cover cap (\*).
- Trendelenburg adjustment up to  $\pm 12^\circ$  (\*).
- Reverse Trendelenburg adjustment up to  $\pm 12^\circ$  (\*).
- High-quality design castors Integral S  $\varnothing$  150 mm, central brake with directional castor.
- Double design pedal to activate the central brake mechanism.

## Mattress platform


- One- or two-section (\*) mattress platform.
- The mattress platform is covered with fire retardant foam, thickness 40 mm and density 45 kg/m<sup>3</sup>.
- The mattress platform is finished with imitation leather.
- The two-section (\*) mattress platform is equipped with an adjustable backrest.
- Adjustable backrest (\*) up to  $\pm 70^\circ$  by means of a gas spring.

## Manipulations

- High/low adjustment by means of pedals that are within easy reach on both sides of the stretcher.
- The Trendelenburg position is activated by means of the foot pedal at the foot end of the stretcher (\*).
- The gas spring-supported backrest can be adjusted by means of a bracket.

## Sideguards Atmosphère (\*)

- Sideguards offer additional protection while moving the patient.
- The Atmosphère sideguards constitute one unit and can be lowered by means of an easy locking mechanism.

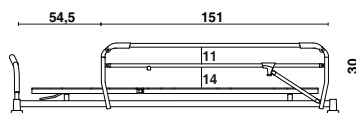

### Sideguards

- Atmosphère (\*)

### Accessories


We would like to refer to our accessorieelst.

(\*) Option



## Dimension

- Width: 79 cm, Length: 202 cm
- Height variation of the mattress platform:
  - 1-section: 43 cm - 83 cm
  - 2-sections: 46 cm - 86 cm
- Safe working load: 250 kg

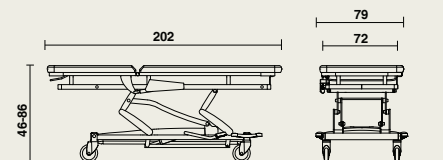

## Materials and finishing

- **Metal:** steel, epoxy coating *RAL 9006 alu*.
- **Synthetic.**
- The hinge points are equipped with self-lubricating aluminium bearings (= maintenance-free).
- Withstand the common cleaning products.

For further information about used materials, constructions and maintenance: consult our materialfiles.

## Options

- Atmosphère sideguards
- (Reverse-) Trendelenburg.
- Fixed oxigen bottle holder.
- Paper roll holder 60 cm.
- Removable tray at head end.
- Removable synthetic cover cap.

[Order here!](#)

**PRAXISDIENST**  
Medical Supplies since 1953

